

USING YOUR TWINJECT[®] AUTO-INJECTOR


STEP 1 - MAKE SURE THAT THE MEDICINE IS READY

Examine your Twinject[®] 0.3 mg or Twinject[®] 0.15 mg regularly. It may not work if the medicine looks cloudy (has particles), pinkish, or more than slightly yellow, or if the expiration date has passed.

In the event of a life-threatening allergic reaction, you should use an out of date product, if that is all you have. Do NOT remove the GREEN caps until you are ready to use your Twinject.

STEP 2 – FIRST DOSE

1. Pull off GREEN end cap marked [1] to see a RED tip. Never put thumb, finger or hand over the RED tip.

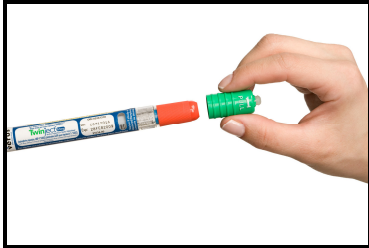

2. Pull of GREEN end cap marked [2].



3. Place RED tip against mid-thigh (injection can go through clothes).



4. Press down hard until the needle enters your thigh through your skin – hold while slowly counting to 10.
5. Remove auto-injector and check RED tip; if needle is exposed, you have received the dose. If not, repeat steps 3 and 4.
6. Prepare for second dose.

Call 911 or the closest emergency department as soon as you administer the first dose.

STEP 3 – PREPARE FOR SECOND DOSE

1. Unscrew and remove RED tip. **Beware of the exposed needle.**
2. Holding BLUE hub at needle base, remove syringe from barrel.



3. Slide YELLOW collar off plunger.



4. **PAUSE HERE.** If your symptoms have not improved within 10 minutes since the first injection, proceed with steps 5 and 6.
5. Put the needle into your thigh (upper leg).



6. Push the plunger down all the way.

Give the used Twinject to a healthcare worker for proper disposal. Do not throw it away in the trash.

Twinject[®]
auto-injector
(epinephrine injection, USP 1:1000)
Available as: 0.3 mg 0.15 mg